




SIS TYPE 30 / 1200 GL AS PER DIN EN 14470-1 (TYPE 30)

Approved for the restricted storage of flammable liquids at the workplace as per DIN EN 14470-1 (type 30) and TRGS 510, annex 1

- ◆ Whole construction tested and approved by an authorized materials testing institute as per DIN EN 14470-1 (type 30) with a fire resistance >30 minutes
- ◆ GS-approved safety (GS-sign)
- ◆ Body made from B1 material RAL 7035 - light grey, doors either in RAL 7035 - light grey or RAL 1007 - daffodil yellow with fire safety glass
- ◆ Height-adjustable shelves at intervals of 32 mm
- ◆ Automatic thermal closing mechanism if fire occurs
 - Door after temperature increase to $\leq 50^{\circ}\text{C}$ and
 - Inlet and outlet channels after temperature increase to 70°C ($\pm 10^{\circ}\text{C}$)
- ◆ DN 80 ventilation connection for exhaust air extraction on cabinet roof
- ◆ Operation complies with TRGS 510 (approved even without exhaust air connection!)
- ◆ Ground connection in accordance with TRBS 2153 (avoidance of sources of ignition)
- ◆ Markings according to DIN EN 14470-1 / TRGS 510 / ISO 3864 / ASR A1.3
- ◆ Leveling feet for installing cabinet horizontally
- ◆ Doors with cylinder lock

With inner fittings:

- ◆ 1 x base sump tray (sump capacity 33 liters) as per WHG - **Ü-sign according to StawaR**
- ◆ 1 x perforated base tray cover as 1st storing level
- ◆ 3 x shelves
- ◆ **door hold-open device**

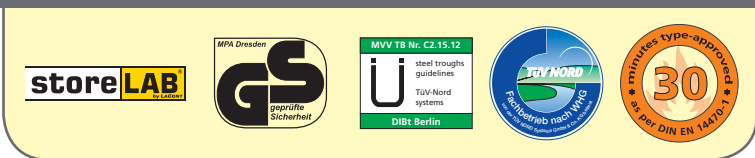

SiS Type 30 / 1200 GL, Article no. B80-2656-A, doors RAL 7035 - light grey, incl. 3 shelves, base sump tray and perforated cover plate



SiS Type 30 / 1200 GL, Article no. B80-2657-A, doors RAL 1007 - daffodil yellow, incl. 3 shelves, base sump tray and perforated cover plate



Socle included, with removable panel



Extraction fans please see page 136!

| Model | | SiS Type 30 / 1200 GL |
|---|--------------------------------|-----------------------|
| External dimensions W x D x H (mm) | | 1195 x 595 x 1935 |
| Internal dimensions W x D x H (mm) | | 1075 x 515 x 1600 |
| Doors with fire safety glass | | 2 |
| Weight (kg) | | 298 |
| Outer body: RAL 7035 - light grey | Material inner fittings | |
| Article no., doors: RAL 7035 - light grey | sheet steel | B80-2656-A |
| Article no., doors: RAL 1007 - daffodil yellow | sheet steel | B80-2657-A |

| Accessories | Material | Article no. |
|-------------|-------------|-------------------|
| Shelf | sheet steel | B80-2610-A |