



Explanations see practical informations on page 6 and 84!


AUS TYPE 30 / 600 AS PER DIN EN 14470-1 (TYPE 30)

Approved for the restricted storage of flammable liquids at the workplace as per DIN EN 14470-1 (type 30) and TRGS 510, annex 1

- ◆ Whole construction tested and approved by an authorized materials testing institute as per DIN EN 14470-1 (type 30) with a fire resistance >30 minutes
- ◆ GS-approved safety (GS-sign)
- ◆ Body made from B1 material RAL 7035 - light grey, door either in RAL 7035 - light grey or RAL 1007 - daffodil yellow
- ◆ Height-adjustable shelves at intervals of 32 mm
- ◆ Automatic thermal closing mechanism if fire occurs
 - Door after temperature increase to ≤50°C and
 - Inlet and outlet channels after temperature increase to 70°C (±10°C)
- ◆ DN 80 ventilation connection for exhaust air extraction in the back wall
- ◆ Operation complies with TRGS 510 (approved even without exhaust air connection!)
- ◆ Ground connection in accordance with TRBS 2153 (avoidance of sources of ignition)
- ◆ Markings according to DIN EN 14470-1 / TRGS 510 / ISO 3864 / ASR A1.3
- ◆ Leveling feet for installing cabinet horizontally
- ◆ Door with cylinder lock

With inner fittings:

- ◆ 1 x base sump tray (BW) as per WHG - **Ü-sign according to StawaR**
- ◆ 1 x perforated base tray cover as 1st storing level
- ◆ **door hold-open device**

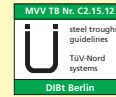

AUS Type 30 / 600 FT, door RAL 7035, Article no. B80-2639-A



Extraction fans please see page 136!



AUS Type 30 / 600 FT, door RAL 7035, Article no. B80-2639-A



AUS Type 30 / 600 FT, door RAL 1007, Article no. B80-2640-A



AUS Type 30 / 600 FT, door RAL 1007, Article no. B80-2640-A

Model		AUS Type 30 / 600 FT
External dimensions W x D x H (mm)		600 x 570 x 600
Internal dimensions W x D x H (mm)		472 x 470 x 425
Sump capacity (BW) (l)		12
Door		1 DIN right
Weight (kg)		63
Outer body: RAL 7035 - light grey	Material inner fittings	
Article no., door: RAL 7035 - light grey	sheet steel	B80-2639-A
Article no., door: RAL 1007 - daffodil yellow	sheet steel	B80-2640-A
Accessories		Article no.
Shelf	sheet steel	B80-2641-A
Socle height 35 mm	B1	B80-2642-A